


Pharmsteri™ PP Cartridge Filters

PSTP series filters provide economical filtration for pharmaceutical and food and beverage applications

PSTP series cartridge filters join our portfolio of filtration solutions that help ensure reliable product processing. Constructed of polypropylene, Pharmsteri cartridge filters are economical and efficient for removing particulate and colloidal contaminants in various pharmaceutical as well as food and beverage applications.

APPLICATIONS

- Pharmaceutical
 - Facility water
 - Buffers
 - LVP/SVP
 - Bulk chemicals
- Food and Beverage
 - Facility water
 - Drinking water
 - Prefiltration for draft beer
 - Clarification for wine/drinks

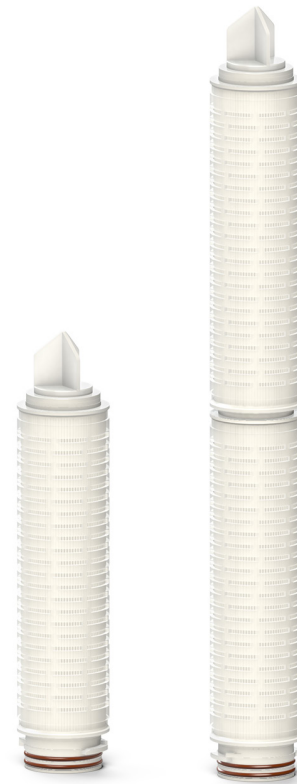

Allows high flow and low pressure, lowering cost of ownership


FEATURES & BENEFITS

All polypropylene construction	Provides excellent chemical compatibility, allowing use in a broad range of applications
Optimized form factor	Allows high flow and low pressure drop, lowering cost of ownership
All components meet USP <88> Class VI Plastics Test	Provides biosafety insurance and global regulatory compliance

SPECIFICATIONS


Retention ratings	0.1, 0.2, 0.45, 1.0, 3.0, 5.0, 10, 20 µm	
Material	Media	Polypropylene
	Support	Polypropylene
	Sleeve, core, top, end caps	Polypropylene
	O-ring	Silicone
	Gasket (Code F)	Silicone
Sealing technology	Thermal bonding without adhesives	
Filtration area	≥0.43 m ² /10" cartridge	
Maximum operation conditions	Maximum forward differential pressure: 4.2 bar @ 25°C (77°F), 2.0 bar @ 80°C (176°F)	
	Maximum operating temperature: 80°C (176°F)	
Multiple sterilization cycles	20 cycles of steam-in-place (SIP) sterilization or autoclave, 30 min at 127°C (261°F)	
Endotoxin releasing	<0.25 EU/mL as determined by the LAL test	
TOC	<0.5 mg/L after a water flush of 11 L/10" cartridge	
Conductivity	≤1.3 µS/cm after water flush of 11 L/10" cartridge	
Particle release	Meets the requirements for WFI set by USP	
NVR	<50 mg/10" cartridge, after autoclave and a 24-hour extraction in 50% ethanol with vibration	
Material toxicity	All components meet the requirement of USP <88> Reactivity Test for Class VI Plastics	
Quality	Manufactured in accordance with ISO 9001 and ISO 13485 certified quality management system	

PERFORMANCE DATA

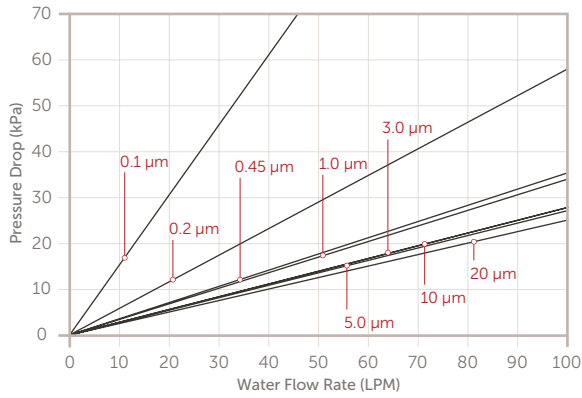


DIMENSIONS


Ø68 mm Code 0 Filter

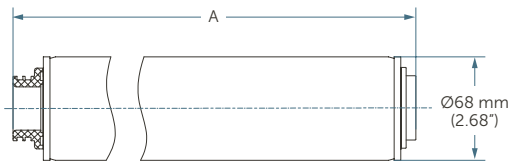


Ø68 mm Code F Filter

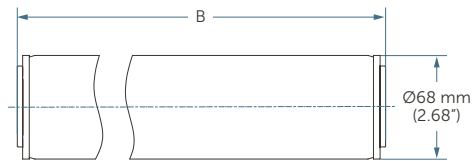

Ø68 mm Code 5 Filter

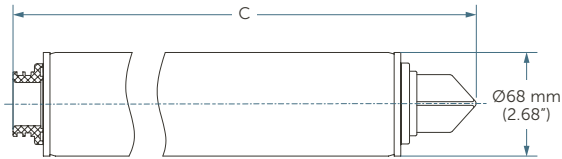


Ø68 mm Code 7 Filter

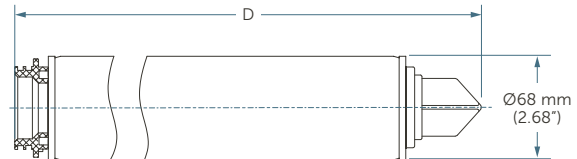

Length	A (Code 0)	B (Code F)	C (Code 5)	D (Code 7)
5"	160 mm (6.3")	137 mm (5.4")	–	–
10"	273.5 mm (10.8")	254 mm (10.0")	314 mm (12.4")	316 mm (12.0")
20"	518.5 mm (20.4")	500 mm (19.7")	559 mm (22.0")	561 mm (22.1")
30"	763.5 mm (30.1")	745 mm (29.3")	804 mm (31.7")	806 mm (31.7")
40"	1008 mm (39.7")	990 mm (39.0")	1048 mm (41.3")	1050 mm (41.3")

ORDERING INFORMATION

Pharmsteri PP Cartridge Filters: part number

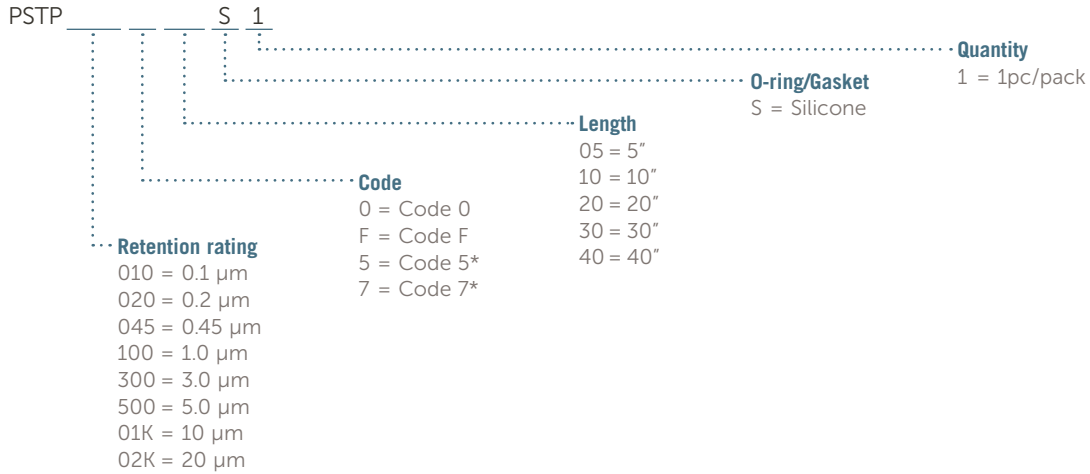

*Code 5 and code 7 are only available with 10", 20", 30", or 40" lengths.

FOR MORE INFORMATION

Please call your Regional Customer Service Center today to learn what Entegris can do for you. Visit [entegris.com](https://www.entegris.com) and select the [Contact Us](#) link to find the customer service center nearest you.

TERMS AND CONDITIONS OF SALE

All purchases are subject to Entegris' Terms and Conditions of Sale. To view and print this information, visit [entegris.com](https://www.entegris.com) and select the [Terms & Conditions](#) link in the footer.



Corporate Headquarters

129 Concord Road
Billerica, MA 01821
USA

Customer Service

Tel +1 952 556 4181
Fax +1 952 556 8022
Toll Free 800 394 4083

Entegris®, the Entegris Rings Design®, and other product names are trademarks of Entegris, Inc. as listed on [entegris.com/trademarks](https://www.entegris.com/trademarks). All third-party product names, logos, and company names are trademarks or registered trademarks of their respective owners. Use of them does not imply any affiliation, sponsorship, or endorsement by the trademark owner.

©2021-2023 Entegris, Inc. | All rights reserved. | Printed in the USA | 9000-11698ENT-0923