

L5P Series Cartridge Filter

Provides flexibility for lab scale to commercial membrane options

L5P cartridge filters are small effective area filters that can be optimized with different media for both liquid and air and gas applications.

The filters can be outfitted with polypropylene (PP) or PES media for small-volume liquid applications. When used with hydrophobic PTFE, these filters become sterilizing-grade and can be used for either air and gas applications or venting.

These single open-end filters have a 56 mm outer dimension, which provides a smaller effective area compared to other 68 mm or 70 mm standard cartridge filters.

APPLICATIONS

- Lab-scale or small volume applications
- Bioburden reduction or sterile filtration (with PES membrane)
- Sterile filtration of air and gas or venting (with PTFE membrane)



Flexibility for both experimentation and production.

FEATURES & BENEFITS

All components meet USP<87><88>, 21 CFR, and commission regulation for food contact materials

Provides biosafety insurance and global regulatory compliance

Multiple media options with different filter lengths

Provides flexibility of utilized in different trials and processes

100% integrity tested before release (filters with PES or PTFE membrane)


Provides safe, sterile filtration for air and gas in pharmaceutical applications

SPECIFICATIONS

Materials of construction	Membrane	PES/PTFE/PP
	Support	PP
	Core and Cage	PP
	End Caps	PP
	O-rings/Gasket	Silicone
	Sealing technology	Thermal bonding, no adhesives
Dimensions	Diameter	Ø56 mm
	Length	2.5", 5"
Filtration area	≥2.2 ft ² per 5" cartridge	
Retention rating	0.1, 0.22, 0.45 µm (Depends on membrane choice, contact Entegris for details)	
Maximum differential pressure	Forward: 4.2 bar @ 23°C (73°F), 1.5 bar @ 85°C (185°F)	
Food contact material	All components meet the requirements of (EC) No 1935/2004	
Indirect food additive	All components meet the FDA Indirect Food Additive requirements cited in 21 CFR 177	
Quality	Manufactured in accordance with ISO 9001 and ISO 13485 certified quality management system	

ORDERING INFORMATION

L5P Series Cartridge Filter: part number

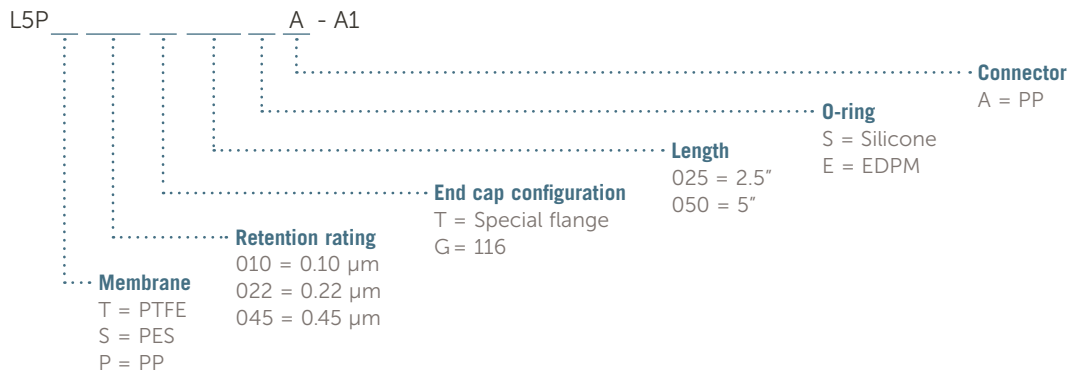

Contact Entegris for details on specific membrane choice, retention rating, and part number.

FOR MORE INFORMATION

Please call your Regional Customer Service Center today to learn what Entegris can do for you. Visit entegris.com and select the [Contact Us](#) link to find the customer service center nearest you.

TERMS AND CONDITIONS OF SALE

All purchases are subject to Entegris' Terms and Conditions of Sale. To view and print this information, visit entegris.com and select the [Terms & Conditions](#) link in the footer.



Corporate Headquarters

129 Concord Road
Billerica, MA 01821
USA

Customer Service

Tel +1 952 556 4181
Fax +1 952 556 8022
Toll Free 800 394 4083

Entegris®, the Entegris Rings Design®, and other product names are trademarks of Entegris, Inc. as listed on entegris.com/trademarks. All third-party product names, logos, and company names are trademarks or registered trademarks of their respective owners. Use of them does not imply any affiliation, sponsorship, or endorsement by the trademark owner.

©2023 Entegris, Inc. | All rights reserved. | Printed in the USA | 4463-13099ENT-0723